



made for you!



JOYCEone

JOYCEone

Your number 1 choice for a perfect fit in no time.

Time-consuming settings are history. Today we offer you JOYCEone, the ventilation mask in one size that fits nearly everyone. This mask, the result of years of experience and mask expertise, is available in a nasal or full face model.

Hand-in-hand: comfort and effective therapy

The forehead support's spring action automatically ensures that the proper distance is maintained for patient comfort and effective therapy. The single layer mask cushion is just as cleverly designed. It finds the ideal position on its own, regardless of the size and shape of the patient's nose and face.

Design is in the details

Your patients will love JOYCEone! This mask features more than a perfect fit and wearing comfort. The unmistakable color-coding of headgear and clips and the convenient recessed grip help your patient to assemble JOYCEone quickly and correctly right from the start.

Modern Efficiency

Short and sweet: the intelligent JOYCEone helps you to save time without compromising on patient care. JOYCEone – made for you!

The advantages of JOYCEone

- Saves time: automatic adjustment of forehead support **3** and mask cushion **1**
- Unmistakable assembly: color-coding on headgear and clips **4**
- Economical: one size fits nearly all patients
- Convenient: effortless cleaning
- Pleasant: quiet exhalation system
- Stable: very durable materials



JOYCEone Nasal

Mask fit on different nose sizes



JOYCEone Full Face

Mask fit on range of nose and face sizes

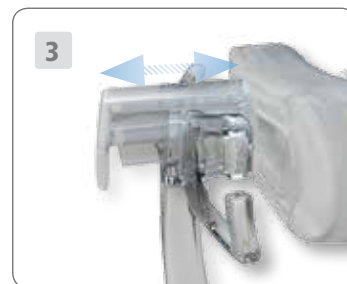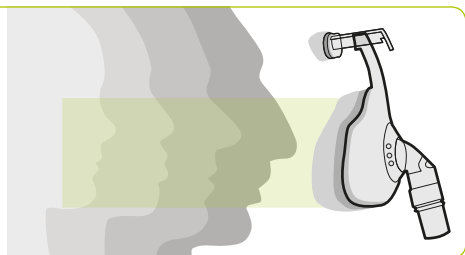

- 1** Mask cushion
- 2** Forehead cushion
- 3** Forehead support
- 4** Headgear clips
- 5** Hook-and-loop tab
- 6** Pressure measurement port
- 7** Quick-release cord (only JOYCEone Full Face)



JOYCEone Nasal mask



JOYCEone Nasal



JOYCEone Full Face
Full face mask

Accessories



Headgear
incl. headgear clips
WM 25255




Headgear clips
2 pieces
WM 25257



Forehead cushion
silicone,
WEINMANN SilkTec
coating
WM 26801



Forehead support
JOYCEone Nasal
incl. silicone forehead cushion,
WEINMANN SilkTec coating
WM 15907
JOYCEone Full Face
incl. silicone forehead cushion,
WEINMANN SilkTec coating
WM 15925

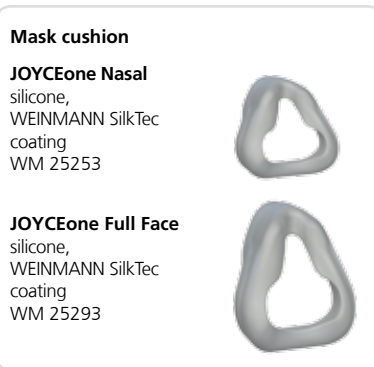

Mask cushion
JOYCEone Nasal
silicone,
WEINMANN SilkTec
coating
WM 25253
JOYCEone Full Face
silicone,
WEINMANN SilkTec
coating
WM 25293



**Connection for pressure
measurement port**
WM 26849



Rotating sleeve
JOYCEone Nasal
WM 26846
JOYCEone Full Face
WM 26847



Ball-and-socket set,
JOYCEone Nasal
retaining ring,
elbow,
rotating sleeve,
plug
WM 15908
JOYCEone Full Face
retaining ring,
elbow,
anti-asphyxia valve,
rotating sleeve,
plug
WM 15944



Carrying bag
WM 26295

Our complete selection of therapy solutions, accessories and mask systems can be found at: weinmann-medical.com

Technical Data


	JOYCEone Nasal	JOYCEone Full Face	JOYCEone Nasal	JOYCEone Full Face
Product class as per Directive 93 / 42 / EWG	II a			
Dimensions (W x H x D)	96 x 115 x 83 mm	120 x 150 x 100 mm		
Weight	87 g	100 g		
Therapy pressure range	4 hPa to 30 hPa			
Tube connection Tapered connection complies with EN ISO 5356-1	Ø 22 mm (male)			
Pressure measurement port	Ø 4 mm			
Temperature range ■ Operation: ■ Storage:	+5 °C to +40 °C -20 °C to +70 °C			
Flow resistance at 50 l/min at 100 l/min			0.05 hPa 0.28 hPa	0.02 hPa 0.28 hPa
Flow resistance Anti-asphyxia valve			50 l/min ≤ 0,5 hPa 100 l/min ≤ 2,0 hPa	
Quoted two-figure noise emission value as per ISO 4871 ■ Sound pressure level: ■ Sound power level: ■ Uncertainty factor:			15 dB(A) 23 dB(A) 3 dB(A)	
Service life			up to 12 months ¹	

¹ Materials used in the manufacture of masks age when exposed to aggressive cleaning products. Therefore, it may sometimes be necessary to replace mask parts earlier.

Product name	WM no.
JOYCEone Nasal	WM 25260
JOYCEone Full Face	WM 25290



83699-EN-07/2014 © Copyright Weinmann, Hamburg. Duplication of any kind only with the express permission of Weinmann. All rights to design and specification modifications reserved.

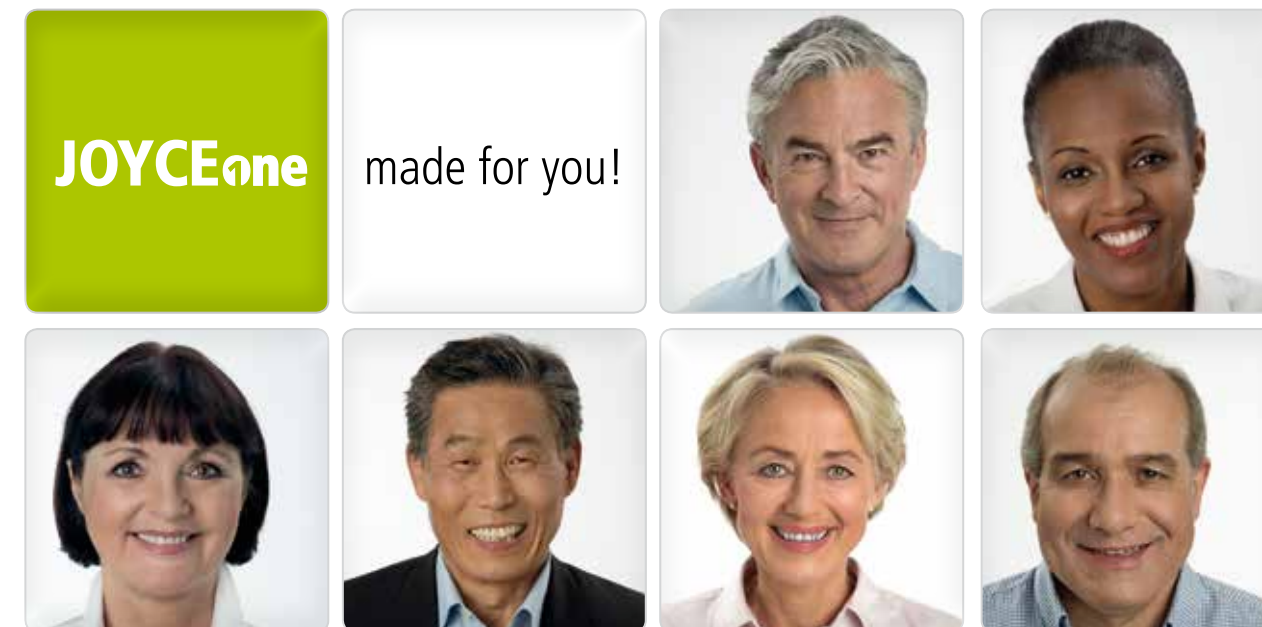

JOYCEone

Full Face